



Delivering the Right Results

Predicting the Initial Deformation and Final Temperatures of Coal Ash Fusibility using the Measured Softening and Hemispherical Temperature

Kim Marshall, Mason Marsh, Trace Yates of LECO Corporation, Saint Joseph, Michigan and John T. Riley, Bowling Green, Kentucky

INTRODUCTION


Ash fusion temperatures give an indication of the softening and melting behavior of fuel ash. Ash fusion temperatures are useful in predicting the slagging and fouling characteristics of fuel ash. The ASTM International standard test method for determining ash fusion temperatures (ASTM D 1857) uses four temperatures to characterize the fusion behavior of coal and coke ashes. The four transition temperatures are as follows:

Initial deformation temperature (IT) Temperature at which the first rounding of the apex of the cone occurs.

Softening temperature (ST) Temperature at which the cone has fused down to a spherical lump in which the height is equal to the width of the base.

Hemispherical temperature (HT) Temperature at which the cone has fused down to a hemispherical lump at which point the height is one-half the width of the base.

Fluid temperature (FT) Temperature at which the fused mass has spread out in a nearly flat layer with a maximum height of 1/16 in.

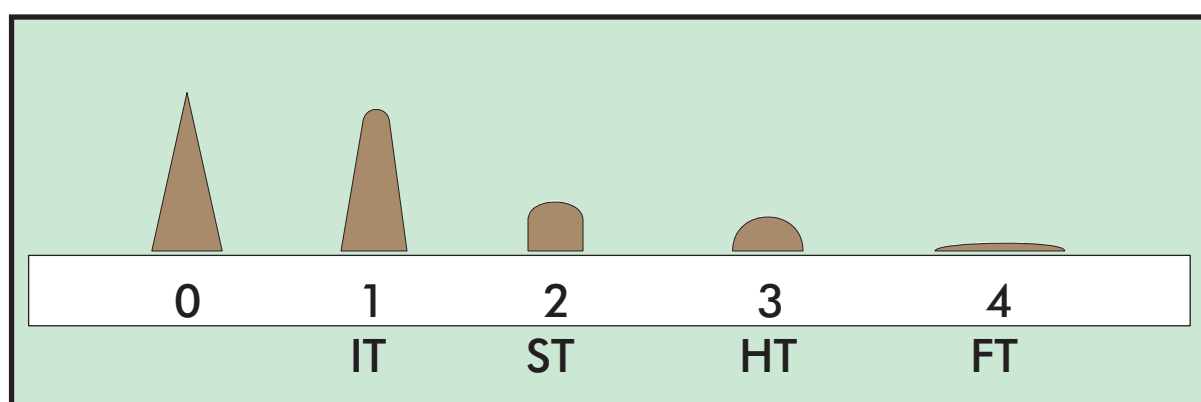

The temperature measurements are made by observing the behavior of 1/4 x 3/4 in. ASTM cones or DIN pins prepared from ash and mounted on a ceramic wafer. The observations may be made under either oxidizing (air) or reducing (60% CO and 40% CO₂) atmosphere. The samples are heated rapidly to 800°C and then, at a rate of 8°C per minute (± 3°), to a maximum temperature (typically 1550°C). The four temperatures are assigned according to the definitions.

The AF700 Ash Fusion Determinator software automatically determines ash fusibility temperatures using Image Recognition Functions (IRF) within the software. Difficulties in determining IT and FT led to the investigation of using the measurable ST and HT values to determine IT and FT. A patent-pending series of equations was developed and verified with ash fusion temperatures from coal and coke ashes from widely varying sources on ASTM cones, DIN pins, and ASTM cones with broken tips.

METHODOLOGY




Equipment Required

- LECO AF700 Ash Fusion Determinator
- ASTM Cone mold or DIN pin mold

Procedure

1. Prepare either ASTM cone or DIN pin for analysis.
2. Place wafer with samples into AF700.
3. Analyze.

Predictive Algorithm Development Procedure

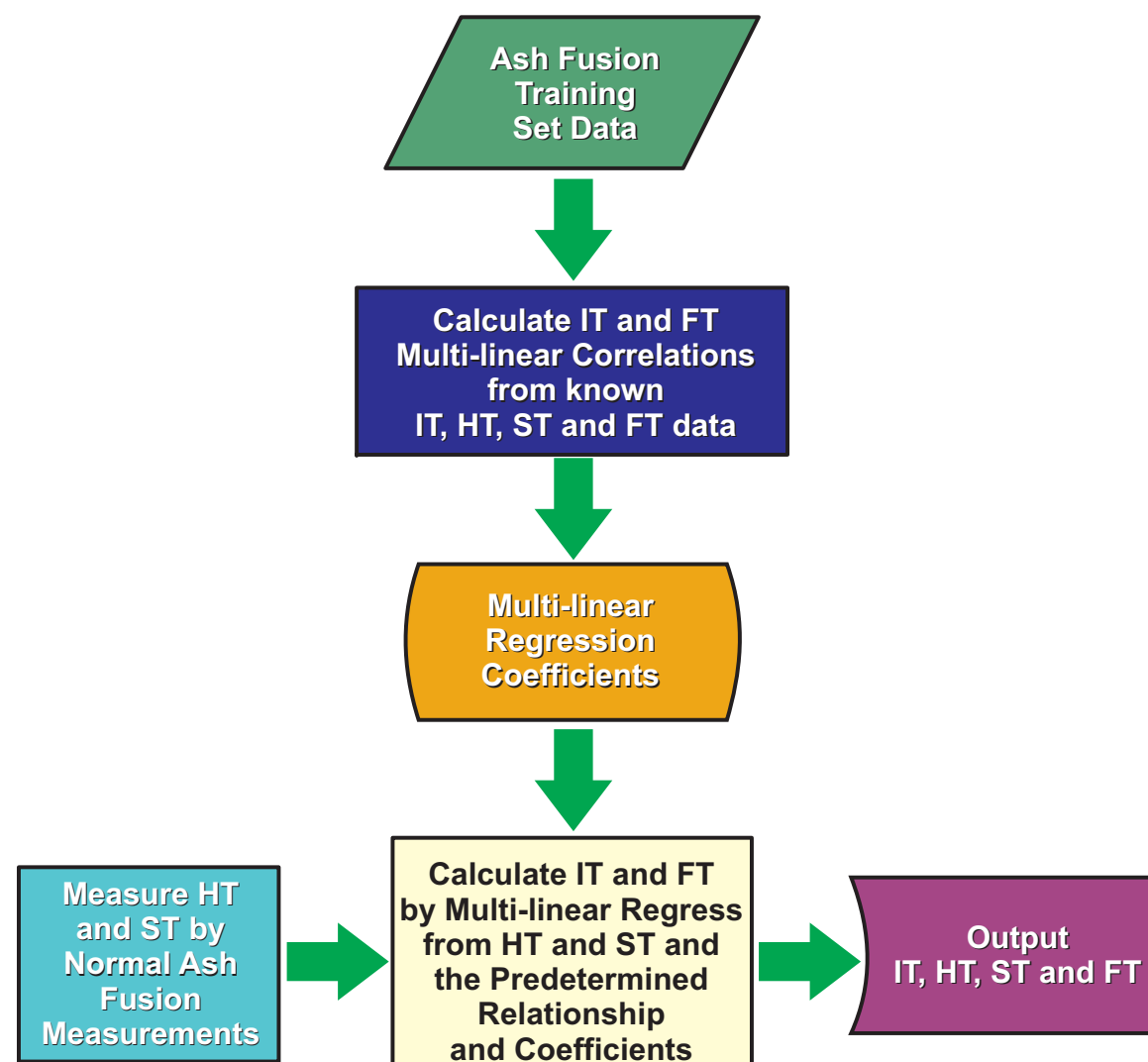

Typical ash fusion determination requires a well defined ASTM cone for accurate IT measurements. (see Figure 1)

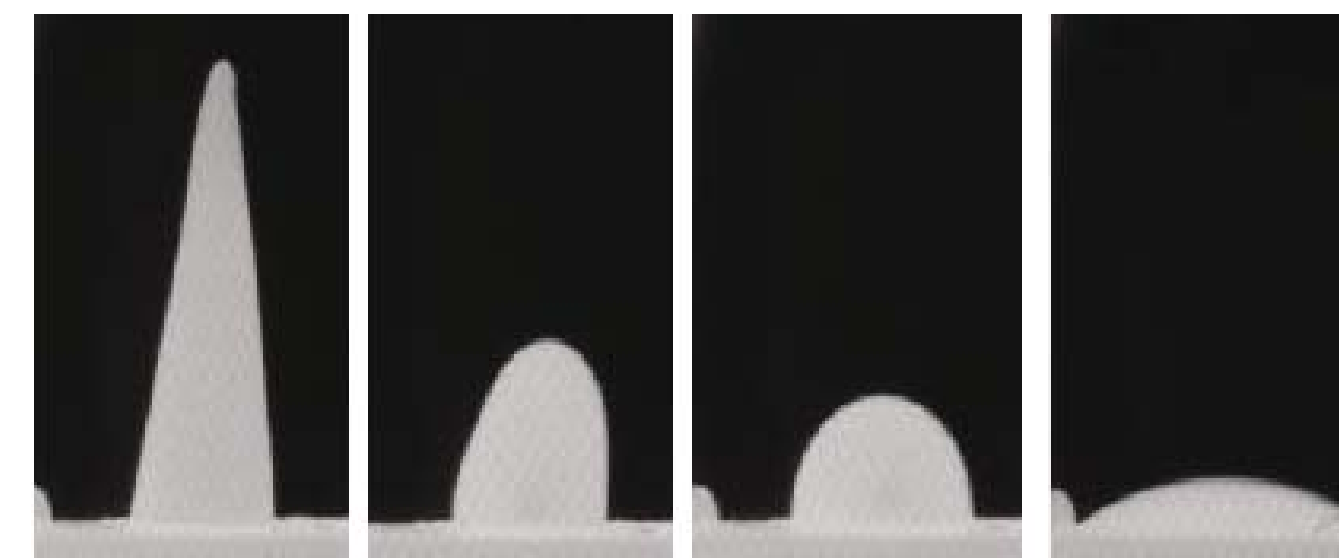


Figure 1. Sample 2000-1 ASTM cone

However, using the patent pending method, both DIN pins (Figure 2) and ASTM cones with broken tips (Figure 3) may be used to determine ash fusion values.

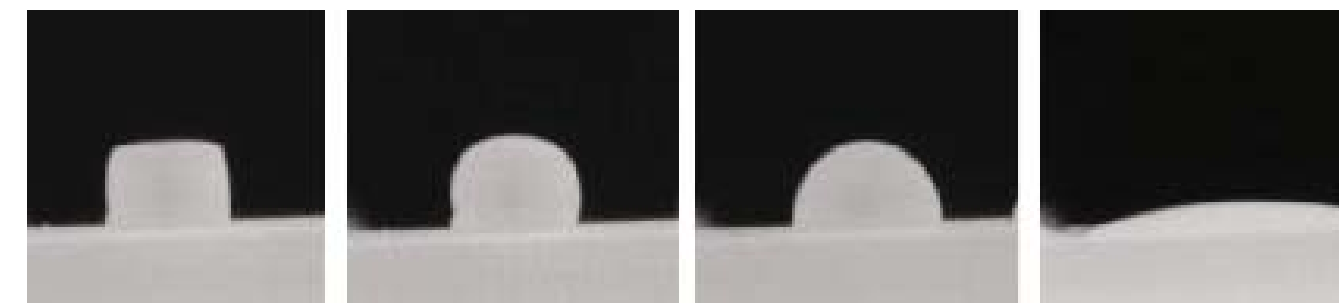


Figure 2. Sample 2000-1 DIN pin

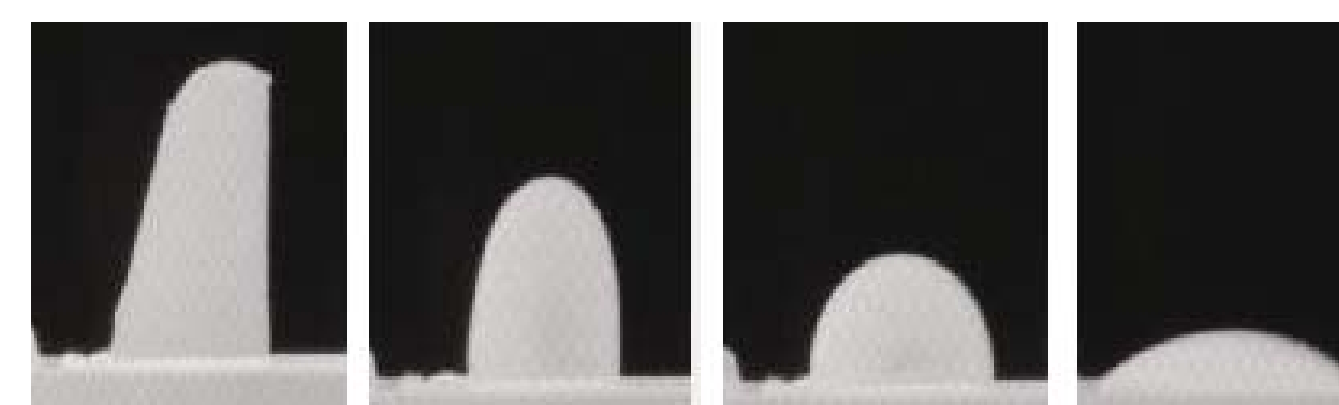


Figure 3. Sample 2000-1 ASTM cone with broken tip

Table 1. Measured vs. Calculated for Various Cone Shapes

Sample	IT ASTM	IT Predicted	ST	HT	FT ASTM	FT Predicted
2000-1 broken cone		1136	1156	1180	1270	1231
2000-1 pin	1072	1129	1151	1179	1284	1235
2000-1 cone	1077	1102	1136	1184	1272	1259
2000-1 pin		1132	1152	1176	1287	1230
2000-1 cone	1072	1127	1150	1179	1263	1236
average	1073	1125	1149	1180	1275	1238
Riley	1104		1147	1175	1225	
Canspex	1096	1117	1138	1165	1223	1222
std	2.4	13.4	7.6	2.9	10.0	11.9
rsd	0.22	1.19	0.66	0.24	0.79	0.96

RESULTS

Table 2. CANSPEX Data With Associated Manually Calculated Predicted Values

Sample	Measurement #	Labs	IT		FT	
			Consensus	Predicted	Consensus	Predicted
2000-1*	16	28	1096	1117	1223	1222
2001-1*	16	26	1104	1104	1193	1195
2001-2*	16	36	1216	1208	1333	1317
2001-3*	16	22	1111	1107	1256	1260
2001-4*	16	24	1215	1213	1253	1269
2002-1*	16	21	1134	1130	1267	1277
2002-2*	16	23	1221	1214	1294	1280
2002-3*	16	24	1247	1236	1340	1322
2002-4*	16	21	1114	1110	1284	1322
2003-1*	16	22	1282	1269	1362	1342
2003-2*	16	20	1330	1315	1411	1391
2003-3*	16	21	1184	1205	1355	1350
2003-4*	16	19	1322	1313	1474	1476
2004-1*	16	21	1433	1432	1559	1553
2004-2*	16	19	1526	1487	1542	1544
2004-3*	16	22	1099	1110	1178	1182
2004-4*	16	22	1183	1176	1420	1376
2005-1*	16	19	1390	1392	1486	1488
2005-2*	16	24	1159	1157	1229	1237
2005-3*	16	22	1106	1109	1199	1200
2005-4*	16	24	1249	1241	1289	1302
2006-1*	16	24	1207	1207	1243	1260
2006-2*	16	22	1152	1158	1334	1324
2006-3*	16	23	1302	1292	1441	1437
2006-4*	16	19	1249	1234	1370	1339
2007-1*	16	29	1134	1139	1292	1294
2007-2*	16	29	1171	1176	1215	1233
2007-3*	16	28	1159	1165	1214	1236
2007-4*	16	27	1390	1378	1493	1496
2008-1*	16	30	1248	1236	1343	1327

*Reducing environment **Oxidizing environment

Table 3. LECO Comparative Study Using AF700 Integrated Prediction Software*

Sample	ST		HT		IT			FT				
	Reference Value (Riley)	AF700	Reference Value (Riley)	AF700	Reference Values (Riley)	CANSPEX	AF700 Predicted	Reference Values (Riley)	CANSPEX	AF700 Predicted		
1997-4	1291	1306	1339	1332	1207		1229	1269	1387	1415	1374	
1998-1	1431	1434	1475	1458	1383		1362	1383	1510	1492	1489	
1998-2	1164	1147	1175	1155	1163		1237	1139	1197	1236	1194	
1998-3	1296	1316	1333	1375	1250		1248	1256	1362	1394	1444	
1998-4	1230	1236	1245	1264	1214		1205	1204	1308	1387	1314	
1999-1	1153	1148	1168	1161	1138		1134	1136	1213	1230	1204	
1999-2	1284	1269	1329	1325	1243		1237	1215	1406	1445	1421	
1999-3	1193	1170	1204	1176	1177		1161	1160	1244	1203	1212	
1999-4	1166	1179	1215	1203	1170		1128	1157	1244	1245	1252	
2000-1	1147	1149	1175	1180	1104	1096	1073	1125	1225	1223	1275	1238
2000-2	1431	1423	1461	1457	1377		1358	1368	1487	1489	1497	
2000-3	1170	1207	1186	1235	1154		1137	1179	1218	1306	1280	
2000-4	1282	1263	1308	1281	1262		1246	1235	1374	1398	1320	
2001-1	1123	1117	1144	1134	1107	1104	1092	1106	1219	1193	1214	1183
2001-4	1224	1230	1234	1234	1208	1215	1222	1215	1255	1253	1251	1263
2002-1	1164	1168	1207	1218	1127	1134	1126	1129	1265	1316	1335	1293
2002-2	1230	1232	1245	1242	1214	1221	1216	1213	1308	1294	1336	1276
2002-3	1265	1264	1284	1302	1245	1247	1244	1223	1340	1340	1383	1358
2002-4	1165	1129	1230	1175	1121	1114	1101	1096	1286	1284	1360	1250
2003-2	1348	1345	1363	1359	1325	1330	1318	1312	1413	1411	1447	1385
2003-4	1373	1371	1415	1403	1327	1322	1335	1321	1469	1474	1471	1446
2004-2							no melt					
2004-3	1121	1115	1138	1134	1105	1099	1072	1102	1189	1178	1207	1186
2004-4	1234	1194	1301	1288	1179	1183	1144	1123	1419	1420	1489	1398
2005-1	1439	1420	1460	1445	1362	1390	1369	1371	1485	1486	1473	1476
2005-3	1133	1116	1155	1148	1109	1106	1102	1095	1221	1199	1224	1211
2005-4	1257	1255	1268	1262	1243	1249	1246	1235	1288	1289	1273	1291
2006-1	1222	1220	1232	1224	1203	1207	1216	1206	1249	1243	1244	1254
2006-2	1220	1139	1262	1210	1180	1152	1119	1091	1337	1334	1334	1300
2006-3	1340	1326	1386	1365	1296	1302	1267	1277	1441	1441	1445	1417

*Reducing environment

OBSERVATIONS and COMMENTS

As stated, the objective of this study was to confirm that the difficult to measure ash fusion values of IT and FT can be accurately predicted from the more easily measurable ST and HT values. ASTM requires a repeatability limit of 30°C between runs. This means that with 95% confidence, the same operator on the same apparatus should be able to have successive determinations that are within 30°C of one another.

From the data shown, all but three of the predicted IT and FT values are well within the 30°C margin of error. The only sample that was consistently outside of the 30°C margin was sample 2004-2 which failed to fall within the specifications for both oxidizing and reducing environments in both studies. Sample 2004-4 fell outside of the specification when analyzed under reducing conditions in the CANSPEX Data evaluation, but was within 30°C of the consensus value when reanalyzed during the LECO Comparative Study. An average absolute difference of 13°C between consensus CANSPEX values and the LECO Comparative Study using the determination values of the AF700 software was observed. The CANSPEX Data with manually determined values showed an absolute average difference of 12°C under reducing conditions and 9°C under oxidizing conditions.

Overall, the predictive algorithm has shown to be a highly accurate tool for determining IT and FT from ST and HT values. The most highly valuable and interesting aspects are the ability of the calculations to account for defects such as broken tips on ASTM cones and the ability to use DIN pins. Figures 1-3 show a standard "good" ASTM cone (Fig. 1), a DIN pin (Fig. 2), and an ASTM cone with a broken tip (Fig. 3). The associated data can be seen in Table 1. Under "normal" conditions, an ASTM cone with a broken tip would be discarded resulting in a loss of operator time, and sample material because IT could not be determined. However, there was no statistical difference for the calculated vs. CANSPEX values when using broken tips. Some operators also prefer to use DIN pins which are easier to make. The IT is very difficult to determine when using DIN pins, but the calculated values from ST and HT are easily determined with good accuracy.

The new predictive algorithm used in conjunction with the AF700 is well within ASTM guidelines for accuracy and precision when determining IT and FT. By automatically terminating the analysis after the last deformation point has been reached for all samples, the algorithm increases throughput while also saving money by extending furnace lifetime. A variety of sample shapes and sizes may also be used, allowing for greater operator freedom to choose either ASTM or DIN pins. Broken ASTM cones may also be used, saving valuable preparation time.

Acknowledgement – Data employed to validate the predictive methods was taken from the CANSPEX Proficiency Test Program operated by Quality Associates International Canada Ltd.